





Astonishing Scanning Experience

Light,

245g featherweight with improved grip for comfort.

Smooth,

Improved speed and optimal scanned images deliver a smoother scanning experience.

Smooth like silk, regardless of PC specs.

Vivid,

Transmits vivid scan data via the new optical system used in the Medit *i*700.

The Brain of the *i*700 with a Splash of Color

Super Accurate, Super Light, Super Smooth.

The digital workflow begins with accurate data. Medit *i*600 precisely acquires data for effortless digital dentistry.

*High accuracy
Full-arch accuracy: 10.9 µm ± 0.98


*Acquires accurate scan data from anterior segment to the occlusal surface



Easy Does It All

Medit *i*600 can be **directly connected to the PC** using just **one cable**, eliminating limitations of wired devices between chairs and clinics for a wider scope of treatment using scanners.



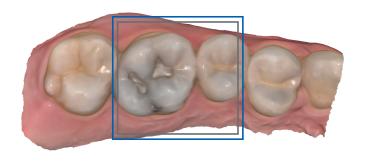


Larger Scan Area (Compared to the i500)

A larger scan area (15X13mm) captures more in one swipe for a more convenient scanning experience.

i600 i500

Making Scanning Easier

Experience amazing digitized data acquisition with **Medit Scan for Clinics**. Our cutting-edge software technology makes scanning more convenient for you.

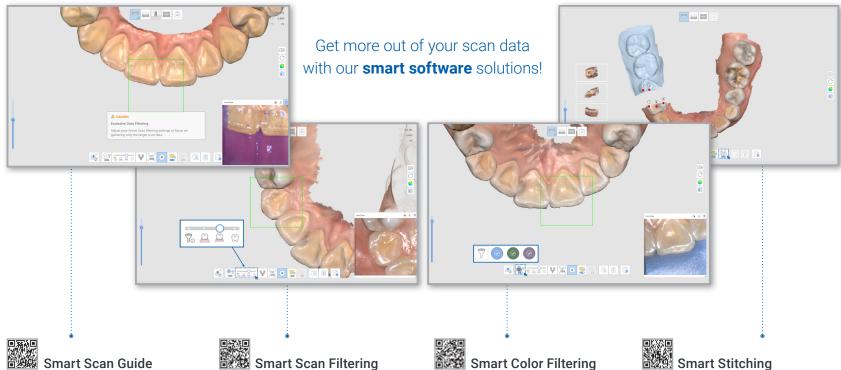

Strengthened user guidance provides help for abnormalities detected during scanning.

Enhanced A.I. filtering performance eases the scanning process by filtering out soft tissue such as lips, cheeks, and tongue.

Automatically filters registered colors, such as the color of medical gloves, during scanning.

Scan with your own scanning strategy. Fragmented scan data are aligned with automatically detected common area.

Increased Patient Satisfaction—

Treatment with accurate data helps improve procedure satisfaction. Consultation using 3D images and Medit's various applications also enhances satisfaction for the overall dental experience.

MEDIT i600

Get More Done on the Chairside

Never worry about the extra cost of using complex programs and systems again! The diversity of Medit Apps makes treatment easier, simpler, and effortless!



into getting an additional scanner, the i600 will make you a perfect partner.

Specifications

| Category | | Medit <i>i</i> 600 |
|---------------------|-------------------------|--|
| Scanning Technology | Scanning Frame | Up to 35 FPS |
| | Imaging Technology | 3D-in-motion video technology 3D full color streaming capture |
| | LED Source | LED |
| | Anti-Fogging Technology | Adaptive anti-fogging |
| Accuracy | Full-Arch | 10.9μm ± 0.98 |
| Handpiece | Dimensions | 248 x 44 x 47.4 mm |
| | Total Weight | 245 g |
| Tip | Tip Size | 22.2 x 15.9 mm |
| | Scan Area | 15 x 13 mm |
| | Autoclavable | Autoclave 100 times 121°C 30 mins 134°C 4 mins (Pre-Vacuum Type) 135°C 10 mins (Gravity Type) |
| | Reversible Tip | Yes |
| Special Mode | Remote Control Mode | No |
| | UV-C Disinfection | No |
| Connectivity | | USB 3.1 Gen1 |

[•] Performance can vary depending on PC specs



[•] Accuracy test of i600 scanner by Medit R&D

[•] Tip for 150 autoclave cycles sold separately

About Us

Medit specializes in 3D measurement and CAD/CAM solutions for dental clinics and labs, including intraoral scanners, based on its own patented state-of-the-art technology. We also develop platform solutions for digital dentistry, supporting collaborative workflows. Additionally, the Solutionix line provides 3D scanners and software to the industrial market. Our goal is to provide innovative technology and the highest quality products to ensure mutual growth for all partners.

Our Medit Community

Continuous feedback, suggestions of new treatment workflows, and active discussion.

From hot debates to fun chats, the Medit community is always lively with our passionate customers.

Colleagues from all over the world actively support and help one another, and sincerely welcome new colleagues into digital dentistry.

Scan the QR code to join!

Medit and its community of colleagues support and wish you the best of luck on your new journey into digital dentistry.

MEDIT i600







User Group

Medit Academy

Help Center



